SERAMON®

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**SUTURE MATERIAL of
POLYTETRAFLUOROETHYLENE**

- **SERAMON® undyed**
- **SERAMON® CHORDAE LOOPS**

DESCRIPTION

SERAMON® is a synthetically manufactured, monofilament, non-absorbable, sterile surgical suture material. The basic suture material consists of polytetrafluoroethylene (PTFE), chemical formula $(C_2F_4)_n$, which has been highly compressed and stretched to achieve a high tensile strength. SERAMON® is only available in undyed form.

The PTFE suture materials are available as combination sutures with or without stainless steel needles and, if required, with specific accessories such as pledgets for cardiac and vascular sutures. They comply with the underlying requirements of the Medical Device Directive 93/42/EEC and apart from marginally deviations with respect to the diameter, the current edition of the European Pharmacopoeia monograph (EP) for "Sutures, sterile non-absorbable (Fila non resorbilia sterilia)".

INDICATIONS

SERAMON® is indicated for soft-tissue adaptation or ligation, in which non-absorbable suture material is indicated, including use in the central circulatory and nervous systems. There has been insufficient experience with use in microsurgery on peripheral vessels or in ophthalmology. PTFE threads can also be used as retaining sutures and for marking purposes.

USE

When selecting and using suture material, the patient's condition, the physician's experience, the surgical technique used and the wound size should be taken into consideration. The established standard knot-tying techniques should be used to ensure knot-holding security. As a result of the even thread surface, additional knots may be necessary in tissues that are under tension. In

order to prevent needle damage during suturing, it is recommended that the needle be grasped at a point between the end of its first third (as seen from its end) and the midpoint of the needle.

ACTION

In wound care, the function of a medical suture is to attach tissue to tissue or tissue to allogeneic and/or xenogeneic material. Ligation or ligature as a special suturing form is used in the closure of hollow organs.

After an initially minimal tissue reaction, suture material made from PTFE is slowly and progressively encapsulated by connective tissue in the body. It is non-absorbable, physiologically inert and thus completely retains its tensile strength. As a result of its monofilament thread structure, SERAMON® can also be used in tissue at risk of infection.

CONTRAINDICATIONS

not known.

WARNINGS / PRECAUTIONS / INTERACTIONS

The demands on the suture material and the techniques required vary according to the therapeutic indication. The user should therefore be familiar with the surgical techniques. When selecting the suture material, its in-vivo properties should be taken into account. As with all foreign bodies, SERAMON® is associated with the risk of calculus formation following long-term contact with the bile ducts or urinary tract. For sutures placed in infected tissues, use of appropriate concomitant measures is recommended. The surgical instruments should be handled carefully, in order to avoid damage to threads. Deformation of needles or holding the needle outside the recommended area can result in needle breakage. Because of the risk of infection from needle pricks, needles should be handled and disposed of with particular care.

ADVERSE REACTIONS

With the correct use of SERAMON®, the following adverse reactions can occur:
transient local irritation, inflammatory reactions

to the foreign body and, in rare cases, formation of suture granulomas or fistulas.

STERILITY

SERAMON® is delivered sterile (sterilization method: ethylene oxide). Do not resterilize, critical changes in suitability for use possible (more details available on request)! Do not use damaged packs! Packs that have been opened for use but then not used are to be discarded!

STORAGE CONDITIONS

No special storage conditions are required. Do not use after the expiry date!

DISPOSAL INFORMATION

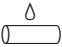
Dispose of unused products or waste material with particular care due to the potential risk of contamination and/or injury. In all other cases observe the specific national regulations.

PRESENTATION

SERAMON® is supplied in a range of gauge sizes and lengths in the form of needle-thread combinations with various needle types or without needles. The combination sutures are available in a sterile pack, individually or with several items in a multipack, as well as furnished with PTFE pledgets to reinforce the body's own structures. More precise details are given in the catalogue.

EXPLANATION OF SYMBOLS ON THE PACKAGE

 1014 CE symbol and identification number of the notified body. The product satisfies the essential requirements of Directive 93/42/EEC on medical devices

 Surgical suture, nonabsorbable, monofilament, undyed

 PTFE Polytetrafluoroethylene

 AR Antireflective needle

 DN Detachable needle

 Pledget (length times width times thickness [mm])

Needle identification symbols and abbreviations are explained in the catalogue.

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SERAG-WIESSNER GmbH & Co. KG